


Horizontal Electrophoresis System (QC-HEP-02Pro)

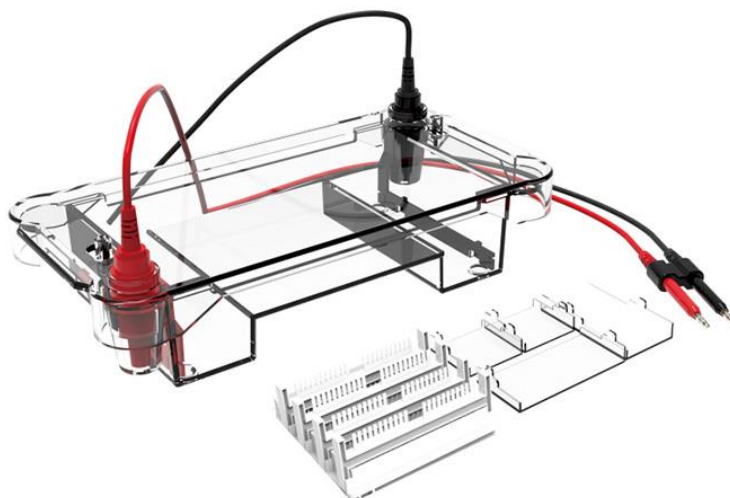

Main Features

1. Polycarbonate is injection molded at one time without leakage.
2. The cooling hole design on the upper cover effectively avoids the condensation of water vapor.
3. Special telescopic elastic glue maker to ensure clean and no leakage of glue.
4. The limit function prevents the positive and negative electrodes from being reversed.
5. Highly flexible wire, power-off design when opening the cover to ensure safety.
6. The electrode head can be replaced for easy maintenance.
7. 99.99% high-purity platinum electrode, the best electrical conductivity, corrosion and aging resistance.

Technical Specifications:

Gel Trays (W*L)	12 cm*12 cm, 12 cm*6 cm, 7 cm*10 cm, 7 cm*7 cm
Buffer Volume	700 ml
Comb Thickness	1 & 1.5 mm
Comb Throughput	2, 3, 6, 8, 11, 13, 18, 25 samples
NO. of Samples	2~50
External Size (LW*H)	270*140*80 mm
Gross Weight	2kg
Standard Accessories	Body tank with electrode*1, Cap*1, Gel tray*5, Gel box*1, 1.5mm thickness 2/3 wells comb*1, 1.0mm thickness 6/13 wells comb*2, 1.0mm thickness 8/18 wells comb*2, 1.0mm thickness 11/25 wells comb*2
Certificate	ISO & CE Certificate